



Fecal Occult Blood Rapid Test




Fecal Occult Blood Rapid Test Cassette (Feces) is a rapid chromatographic immunoassay for the qualitative detection of human occult blood in feces by professional laboratories or physician's offices. It is useful to detect bleeding caused by a number of gastrointestinal disorders, e.g., diverticulitis, colitis, polyps, and colorectal cancer. Fecal Occult Blood Rapid Test Cassette (Feces) is recommended for use in 1) Routine physical examinations, 2) Hospital monitoring for bleeding in patients, and 3) Screening for colorectal cancer or gastrointestinal bleeding from any source.


FEATURES

- Health Canada Class 2 Licensed (106813)
- Rapid results: Results in 5 minutes
- Easy and clear interpretation of results
- Simple to perform sample collection
- No dietary restrictions required
- High sensitivity: Detects as low as 50 ng/mL of hemoglobin or 6 µg/g feces.
- 98% accuracy
- Dual antibody sandwich assay: For accurate and selective detection
- Internal control line: Confirms valid test procedure and specimen adequacy

Specimen Collection:

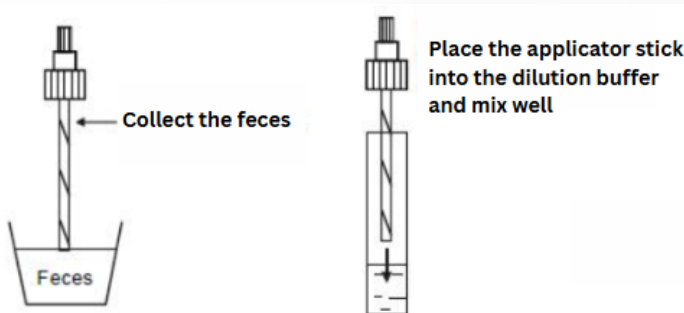

Test Procedure:

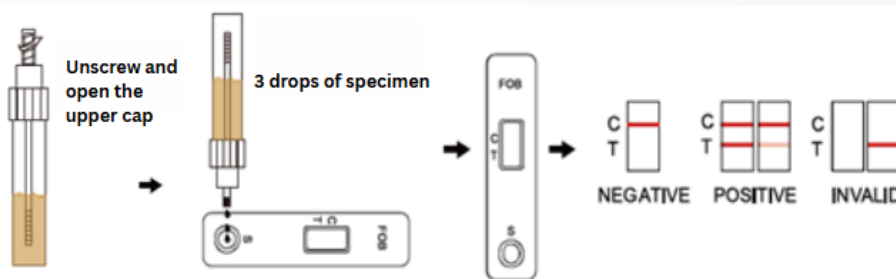

Product Description	Catalog No.	Format	Cut-off Value (ng/ml)	Kit Size
Fecal Occult Blood Rapid Test Cassette (Feces)	GEFOB-602b	Cassette	50	20 Tests/Kit
	GEFOB-602h	Cassette	150	20 Tests/Kit